

### QUADRICEPS TABLE with Backrest

LCS 311

- Quadriceps Exercise chair is very effective and used as a tool in the field of Physiotherapy & Rehabilitation.
- It enables exercise in sitting position i.e. instance knee stretching and exercise in reclining position i.e. knee flexion & hip stretching.
- It has adjustable backrest position, which include a full range of lower and extremity exercises. Backrest can be adjustable with any angle of flexion (0 to 90 Degree).
- Resistance exercise can be adjusted both by increasing & decreasing the weight.



### SHOULDER CUM WRIST EXERCISER

LCS 321




Wall mounting applicable for shoulder & arm exercises. It can be adjusted according to height and resistance can be increased by tightening the knobe, it is made of MS Steel with chrome plated.

### SHOULDER WHEEL

LCS 322

Wall mounting made of MS Steel 100cms wheel dia with calibrated sensitive resistance mechanism. The 360-degree scale enables degree of revolution to be read from either direction. The arc of the attachment can be adjusted at desired height for each patient.

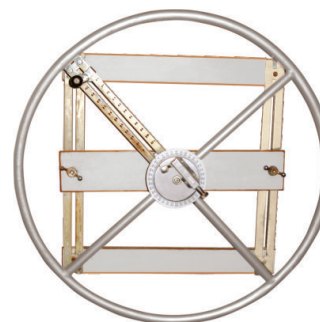

## SHOULDER PULLEY

Made in T shape, suitable for Shoulder & Elbow exercises with two pulleys, rope and handle set.



LCS 331

## SHOULDER ABDUCTION LADDER



Made of teak wood duly polished having 30 numbered steps.

LCS 341

## SUPINATOR PRONATOR EXERCISER


Wall mounted, the circular cylinder frame with grip handle is fitted on the height adjustable bar. The resistance can be increased or decreased with the help of key.



LCS 351

## GRIP EXERCISER

Grip exerciser suitable for weak fingers & hands. This unit is fitted on a laminated wooden base with 6 detachable springs.

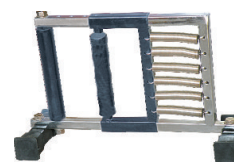

LCS 371

## PEG BOARD ROUND


Wooden laminated board with round slots for peg insertion. Pegs are in 5 different colours suitable for severely handicapped child.



LCS 401

## PEG BOARD SQUARE

Wooden laminated board with square slots for insertion of square pegs. The pegs are in 5 different colours.

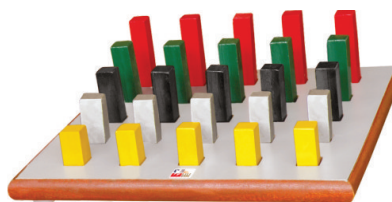

LCS 402

## ANKLE EXERCISER

The foot is mounted on a wooden base with two adjusting springs for Ankle exercises. Tightening spring's nuts can increase resistance further.



LCS 421

## HEEL EXERCISER

LCS 431



The heel unit fitted on a hard laminated wooden base and patient can exercise with sitting on a chair. The pressing rod marked and fitted with a marker.

## STATIC CYCLE

Built in steel tubular frame mounted on plane rubber tipped tubular bar, provided with adjustable seat and specially designed handle available.



LCS 451

## STATIC CYCLE AIR BIKE

500 mm Flywheel for complete exercise.

Permanent Link Dual action handlebars cycling and rowing simultaneously delivering the benefits of exercise to lower body and upper body.

Assembly Size : 1050 x 540 x 1160mm.

Air Resistance : Provide low friction and less wear and tear

Net Weight : 23 Kgs

Max. User Weight : 100 Kgs.



LCS 451H

**CUFF WEIGHT SET:** (A-0.5Kg, B-1Kg, C-2Kg, D-3Kg, E-4Kg & F-5Kg.) each one

LCS 541

Exercise for weak muscles from different weight & sizes.



**DUMBBELL SET:** 1 to 5 kg.

**DUMBBELL:** 1kg.

**DUMBBELL:** 2kg.

**DUMBBELL:** 3kg.

**DUMBBELL:** 4kg.

**DUMBBELL:** 5kg.



LCS 550

LCS 551

LCS 552

LCS 553

LCS 554

LCS 555

Made of Vinyl Dumbbells are available in different sizes and colours 1kg. to 5kg.

**GONIOMETER:** for Fingers


LCS 561

A simple & accurate device for measuring range of motion of small finger joints (I.P. and M.P. joints) constructed of stainless steel. Unit measuring about 14cm. in length



**GONIOMETER:** 180 Degree

LCS 562

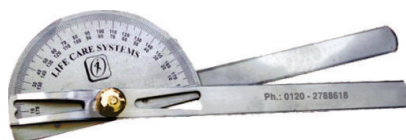


A simple & accurate device to measure joint's range of motion. Stainless steel construction. Two 180° scale in opposite direction for convenient reading of all angles. Unit measures about 50mm x 200mm. in length.

**GONIOMETER:** 360 Degree

LCS 563

Two 360 degree scale projecting in opposite directions around the dial for convenient reading. Aluminium sheet construction with a tension arm knob to hold when desired.

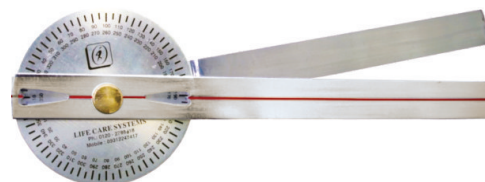

## INFRA RED LAMP: Table Model (150 W)



LCS 571


## INFRA RED LAMP: Floor Model



LCS 572

Luminous bulb (150 W) with portable stand along with dimmer

## TILT TABLE MANUAL

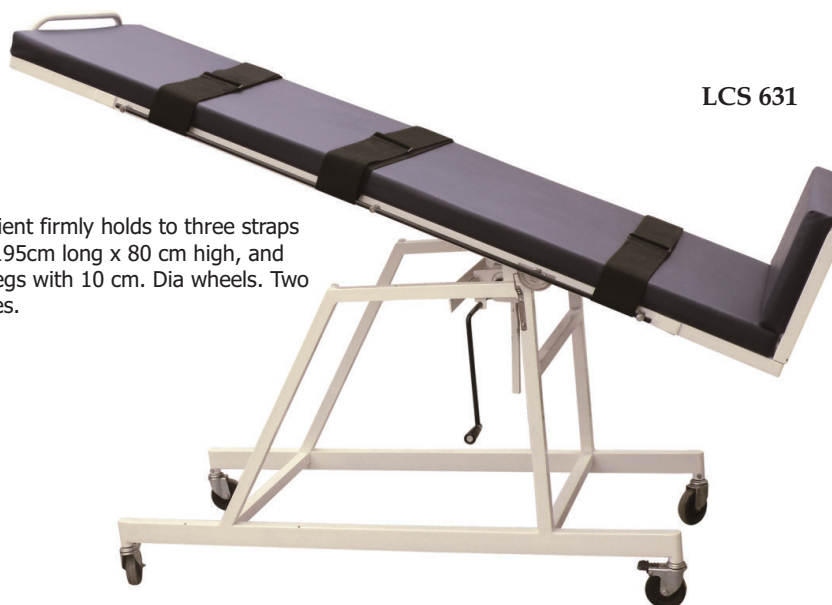

LCS 631

Tilt angles 0 to 90 degree, the patient firmly holds to three straps provided. The tabletop is 61cm x 195cm long x 80 cm high, and on heavy-duty 35 mm steel tube legs with 10 cm. Dia wheels. Two gripping handles for various activities.

LCS 651  
PHYSIOBALL: 55cm.

LCS 652  
PHYSIOBALL: 65cm.

LCS 653  
PHYSIOBALL: 75cm.

LCS 654  
PHYSIOBALL: 85cm.

LCS 655  
PHYSIOBALL: 100cm.

LCS 656  
PHYSIOBALL: 120cm.

